




Quick Guide to TESTING

- 1 Turn the unit ON
- 2 Visually inspect the ear and select an ear tip
- 3 Place the probe in the ear so it fits snugly. Clip the probe to clothing.
- 4 If you would like to use the Child Mode, be sure the icon is selected.

MADSEN Alpha (OAE)

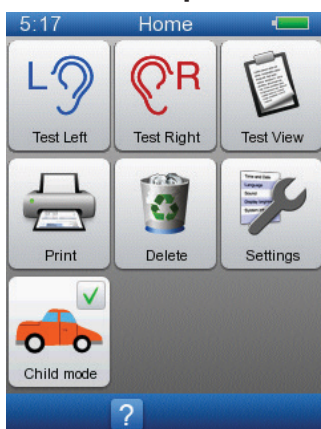


Select the Test Ear.

The speaker check and calibration will begin.

If probe fit is acceptable, the test will automatically begin.

MADSEN Alpha (OAE+)

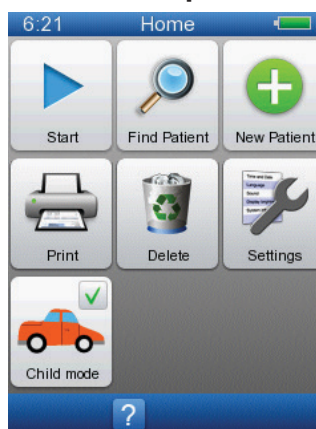


To begin, Select Start or Find Patient (already entered patient) or New Patient (enter a new patient).

Select the desired test protocol and test ear.

The speaker check and calibration will begin.

If the probe fit is acceptable, the test will automatically begin.

At the end of the test, select one of the following

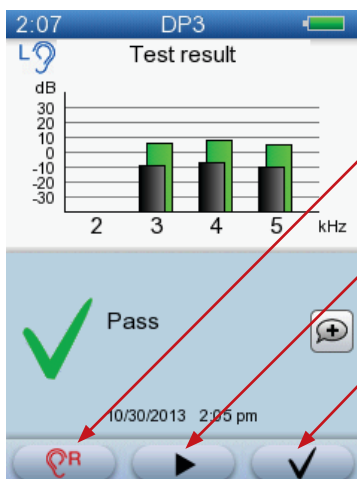

Test the other ear



Retest the same ear

Save/Finished

Quick Guide to the Probe Test

The Probe Test can be started 2 ways:

- 1 Place the probe in the test cavity on the back of the MADSEN Alpha
- 2 From the Main Screen go to Settings, Probe Test and then insert the probe into the test cavity

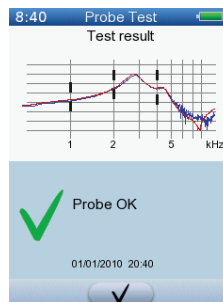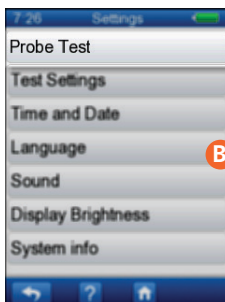

The probe should be placed firmly and straight into the test cavity. If not done properly an incorrect "Fail" reading will be obtained. If this happens, reposition the probe and restart the test.

The speaker check and calibration will begin.