



Elite by Design

Embla Sandman Elite sleep diagnostic software has maintained its “Elite” reputation in the field of sleep diagnostics for over 15 years. This well established and widely accepted software platform delivers the ideal combination of ease and flexibility that is valued by thousands of sleep professionals throughout the world. Sandman Elite is designed to meet the needs of clinical and research facilities large or small.

High Quality

- Exceptional recording, scoring and data management features
- Records high-quality data to produce first-class studies and reports for your laboratory
- Microsoft® Windows 7, 64-bit compatible

Tools You Need

Choose from a complete range of system options to customize tools specific to your laboratory's requirements.

- Select from a range of modules available for computer assisted scoring, data export, and database management
- Bi-directional support of HL7 integration using our exclusive Enterprise™ Practice Management Software
- Craft reports tools with customized formulas, layout, and logo

Intuitive and easy to use Sandman Elite Software is designed with user configurable tools that are available in both Collection and Analysis. Sandman Elite Software offers Real Time On-Screen Digital Video that can be adjusted for increased frame rates during video capture and, when used with the Embla MDrive, will support IP Video collection.

Customizable Toolbars allow for selection of features and tools to streamline workflow and increase collection and scoring efficiency.


Customizable Signal windows and workspaces allow users to view multiple signal windows at different time bases and allow for the simultaneous review of data from different points in time. This useful tool can be used during collection and allows review of look back data and real time data simultaneously.

Digital Display and **Running Digital Displays** of alpha numeric data for any channel or event tag can be customized to view values of signals as they are collected.

Trend Windows and **Event Windows** are customizable and allow for quick and easy navigation to specific areas of the recording. Scoring during Collection or Analysis is performed with the **Single Click Event** placement option that uses channel specific event definitions so events can be added with a single mouse click.

Patient Observations can be customized and allows for collection of observations during acquisition at timed intervals or on demand and is available for review when the study is opened in Analysis.

RDI Summary tool provides information to the recording technologist about the effectiveness of changes in therapy during collection.

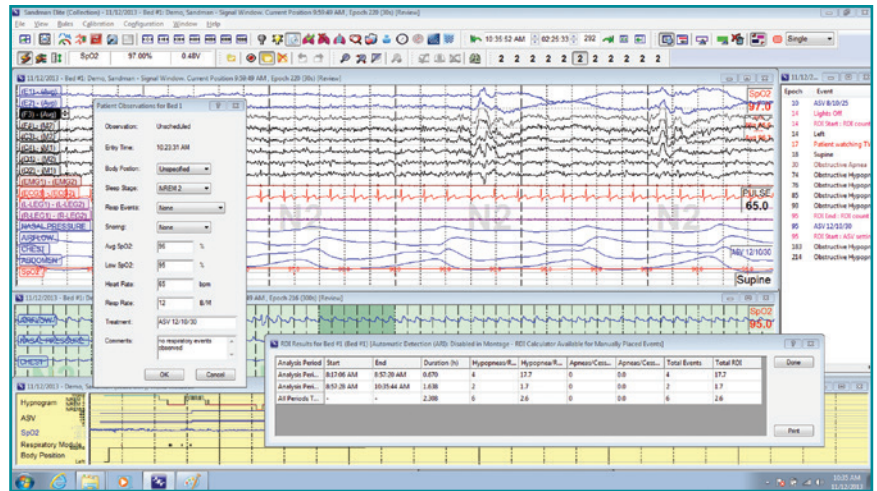

Sandman Elite Recording Tools – Record superior sleep studies

The Embla® Sensor Advantage

- Complete line of sensors and supplies available 24/7 at www.natusneurostore.com
- Unparalleled accuracy with Embla Universal XactTrace® respiratory effort belts
- Optimal sensor cable lengths available

Compatible Amplifiers:

- N7000
- S4500
- SD32 Global
- Embletta Gold



Service

Excellence in customer and technical service is a core value of Natus Neurology.

Essential Support Elements:

- Accessible and effective Technical Support
- Definitive technical documentation and knowledgeable installation teams
- Excellent replacement unit or spare part availability
- Extended warranty and service coverage programs
- Comprehensive, flexible customer training courses

Supplies

Natus Neurology offers a full range of neurodiagnostic supplies, featuring:

- Dedicated and knowledgeable customer service
- Streamlined order processing to save you time
- Convenient online portal (US Customers only)
Natus Neuro Store – www.natusneurostore.com

To learn more about Natus Neurology Service Programs or our full line of Supplies, contact your local distributor or sales representative.

US Customers Call: **1-800-356-0007**

International Customers Call: **+1-608-829-8500**



Natus Neurology Incorporated
3150 Pleasant View Road
Middleton, WI 53562 USA
Tel: 1-800-356-0007
1-608-829-8500
Fax: 1-608-829-8709
www.natus.com

Embla Systems
1 Hines Road, Suite 202
Kanata, Ontario
Canada K2K 3C7
Tel: 1-888-662-7632
1-716-551-6900
Fax: 1-775-860-5726