



Reliable sound for your fitting modules

RoomTune™

Demonstrating hearing instruments is an important part of the fitting process. This is when you set expectations, counsel clients and show what their new hearing instrument can do for them. However, there are many variables to contend with when demonstrating hearing instruments, including the acoustical environment of the fitting room. This one variable can be difficult to control – until now with RoomTune in Aurical Aud.




Most fitting software modules play back the sounds using a dedicated amplifier or speaker that only adjusts the overall level of the sound. RoomTune is new functionality in Aurical Aud that measures how a known sound presents in a given room, and applies an equalizer/frequency compensation during the subsequent playback of sounds. RoomTune provides a more reliable sound for your fitting system so you can give your clients a more realistic and precise demonstration of the hearing instruments.

Better demonstration of hearing instruments

Use Aurical Aud as a sound card

RoomTune functionality is available to any software outside Otosuite® – including fitting modules - as Aurical Aud becomes selectable as a Windows® sound card. Now you can use Aurical Aud and the connected sound field speakers as the compensated sound system with your fitting modules. And, no change to the fitting modules is necessary.

RoomTune is available with Otosuite version 4.75 and higher, and requires the updated RoomTune/Talk back microphone, which is available as an accessory.

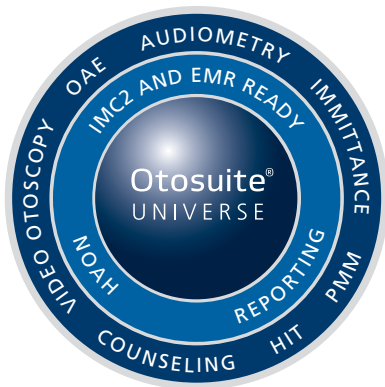

Learn more about our complete fitting solution. Contact your local Natus representative today or visit hearing-balance.natus.com/aurical

KEY BENEFITS

- Easy to use – the RoomTune measurement procedure is automated so there's no need for special training or instructions. Just follow the simple on-screen instructions.
- No need to update fitting modules – use RoomTune with any software outside of Otosuite – Aurical Aud and your connected sound field speakers work as the compensated sound system with your fitting modules.
- More reliable hearing instrument demonstrations

A complete, modular fitting solution, Aurical encompasses audiometry, PMM, HIT, video otoscopy and a HI programming interface in one system. Hearing care professionals can use the modules independently or combine them as needed. This latest release demonstrates once again how Aurical continues to innovate – enabling hearing care professionals to realize greater efficiency gains while boosting the quality of counseling and care.



Healthcare solutions with one thing in mind. You.

©2020 Natus Medical Incorporated. All Rights Reserved. All product names appearing on this document are trademarks or registered trademarks owned, licensed to, promoted or distributed by Natus Medical Incorporated, its subsidiaries or affiliates. 7-26-1640-EN Rev02

natus®

Natus Medical Incorporated

natus.com