



RoomTune. RELIABLE SOUND FOR YOUR FITTING MODULES

Demonstrating hearing instruments is an important part of the fitting process. This is when you set expectations, counsel clients and show what their new hearing instrument can do for them. However, there are many variables to contend with when demonstrating hearing instruments, including the acoustical environment of the fitting room. This one variable can be difficult to control – until now with RoomTune™ in AURICAL® Aud.


AURICAL®


otometrics

BETTER DEMONSTRATION OF HEARING INSTRUMENTS

Introducing RoomTune™

Most fitting software modules play back the sounds using a dedicated amplifier or speaker that only adjusts the overall level of the sound. RoomTune™ is new functionality in AURICAL® Aud that measures how a known sound presents in a given room, and applies an equalizer/frequency compensation during the subsequent playback of sounds. RoomTune provides a more reliable sound for your fitting system so you can give your clients a more realistic and precise demonstration of the hearing instruments.

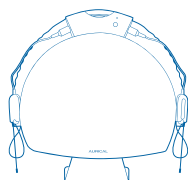

AURICAL – YOUR WAY TO FITTING EXCELLENCE

Use AURICAL Aud as a sound card

RoomTune functionality is available to any software outside OTOsuite® – including fitting modules - as AURICAL Aud becomes selectable as a Windows® sound card. Now you can use AURICAL Aud and the connected sound field speakers as the compensated sound system with your fitting modules. And, no change to the fitting modules is necessary.

RoomTune is available with OTOsuite version 4.75 and requires the updated RoomTune/Talk back microphone, which is available as an accessory.



See more about our complete fitting solution at www.otometrics.com/aurical

KEY BENEFITS

- **Easy to use** – the RoomTune measurement procedure is automated so there's no need for special training or instructions. Just follow the simple on-screen instructions.
- **No need to update fitting modules** – use RoomTune with any software outside of OTOsuite – AURICAL Aud and your connected sound field speakers work as the compensated sound system with your fitting modules.
- **More reliable hearing instrument demonstrations**

Streamline your fitting process with AURICAL®

AURICAL is a complete fitting solution that encompasses audiometry, PMM, HIT, video otoscopy and a HI programming interface in one modular system. You can use the modules independently or combine them as needed. AURICAL can help you to realize greater efficiency gains while boosting the quality of counseling and care in your clinic.

